



FULL SIZE PORTABLE SCOOTER

- Seat Color
Black Vinyl
- Body Colors
Blue
Red

Operating range up to 12 miles!

LiteRider GOLDEN

- 3 wheel GL111D
- 4 wheel GL141D

Enjoy the best of both worlds with the LiteRider™ scooter! The LiteRider is a full-size scooter with all of the benefits of a compact, portable scooter! The LiteRider comes standard with a Delta Tiller, a large 17"x16" stadium style seat, and adjustable armrests. With a 300 lb. weight capacity, the LiteRider can be quickly and wirelessly disassembled for transport or storage.



Model	GL111D	GL141D
Medicare Code	K0800	K0800
Weight Capacity	300 lbs	300 lbs
Drive Wheels	Rear	Rear
Maximum Speed	5 mph	5 mph
Operating Range up to ¹	12 miles	12.5 miles
Ground Clearance Under Deck	3.75"	3.75"
Ground Clearance to Transaxle/Motor	2.5"	2.5"
Turning Radius	39"	48"
Type Batteries	2-22AH	2-22AH
Colors	Red, Blue	Red, Blue
Freewheel Mode	Yes	Yes
Electronic Speed Control	Yes	Yes
Electro-Mechanical Brakes	Yes	Yes
Charger	Off Board	Off Board
	90 Amp	90 Amp
Controller	Dynamic R Series	Dynamic R Series

Dimensions		
Length	41.5"	43.25"
Width	22"	22"
Height (ground to top of back of seat)	37" - 39"	37" - 39"
Ground to Top of Seat	23.5" - 25.5"	23.5" - 25.5"
Ground to Top of Deck	6"	6"
Deck to Top of Seat	18" - 20"	18" - 20"
Number of Seat Height Positions	3	3
Size of Increments	1"	1"
Front Axle to Rear Axle	30"	30"

Weight of Unit: (assembled)	125 lbs	133 lbs
Front Half	44 lbs	52 lbs
Rear Half	28 lbs	28 lbs
Seat w/Arms	21 lbs	21 lbs
Batteries	29 lbs/pack	29 lbs/pack

Tires:	Solid/Low Profile	Solid/Low Profile
Front	9"	9"
Rear	9"	9"
Rear Anti-Tippers	2"	2"

Standard Seat	Stadium Style	Stadium Style
Back Height - no headrest	16"	16"
Width x Depth	17" x 16"	17" x 16"
Color	Black Vinyl	Black Vinyl

Optional Seat	Stadium Style	Stadium Style
Back Height - no headrest	16"	16"
Width x Depth	20" x 18"	20" x 18"
Color	Black Vinyl	Black Vinyl

Warranty:		
Frame	Lifetime	Lifetime
Drive Train	2 Years	2 Years
Electronics ²	2 Years	2 Years
Batteries	1 Year	1 Year

¹ Battery range at 200 lbs. and will vary due to rider weight, drive surface, terrain, and battery type.
² Electronics warranty excludes batteries.

Features


- Delta Tiller standard
- 300 lb. weight capacity
- Rugged anti-rust steel frame tubing
- Rugged plastic floor deck for optimum durability
- Disassembles quickly and easily
- LED ultra bright adjustable headlight



Red



Blue

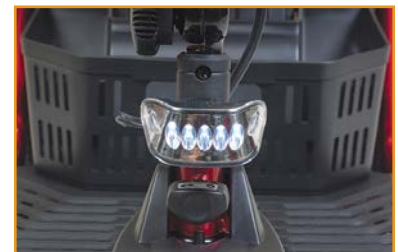

LED ultra bright adjustable headlight

Enjoy quick and wireless disassembly!



Pick up the seat



Pull out the battery pack



Remove the baskets from under the seat and from the tiller (not shown)



To separate the rear transaxle from the frame, squeeze the release lever located on the underside of the frame



While squeezing the release lever, pull up on the rear frame



Disassembles into six pieces for transport and storage! No wires to connect or disconnect!

Literature is current at time of printing. Golden Technologies reserves the right to make changes to the product or literature at any time.